


Acarbose

Method:

Column	Ultisil XB-NH2, 5µm, 4.6×250mm (P/N 00204-31043)
Mobile Phase	Phosphate Buffer:Acetonitrile =25:75
Detection	210nm
Temperature	35°C
Flow Rate	2.0ml/min
Injection Volume	10ul
Note	/

Chromatogram and Data:

System Suitability

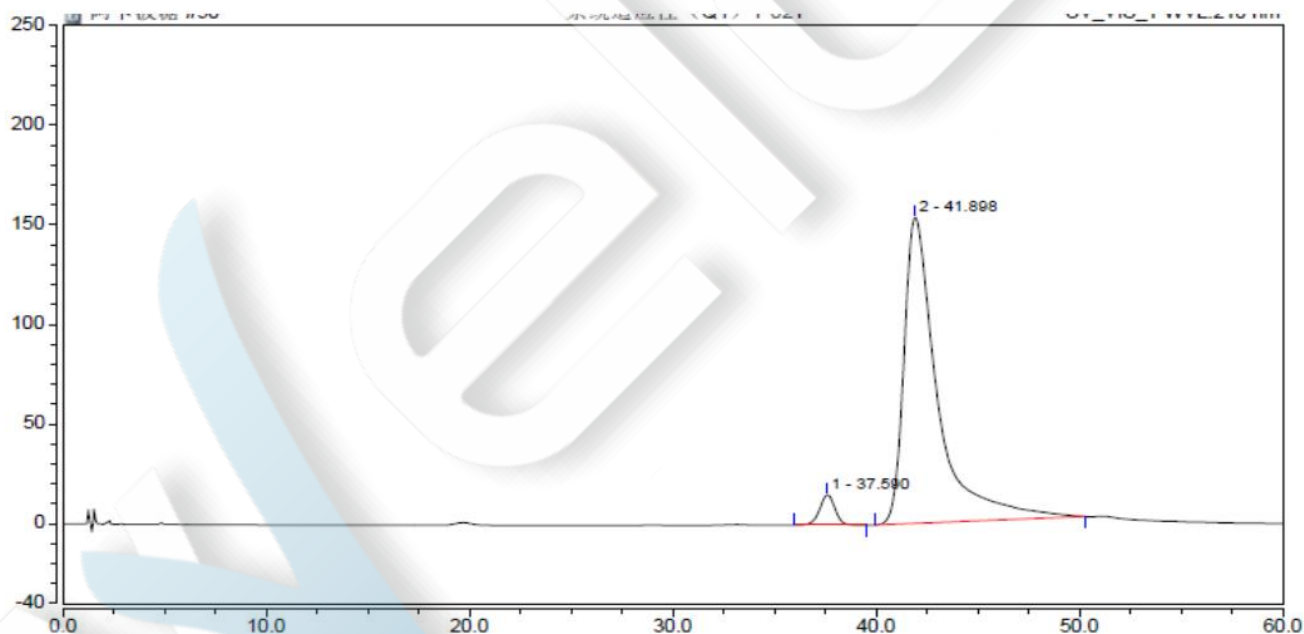

No.	Peak name	Rt min	Area mAU*min	Height mAU	Resolution (EP)	Peak-to-Valley Ratio	Plates (EP)
1		37.590	13.191	15.202	2.21	9602.01	12455
2		41.898	295.966	153.660	n.a.	293989.16	4262

Unless otherwise stated, the results shown in this test report refer only to this sample tested. The report can not be copied without the permission of Welch Materials, Inc.