


Diniconazole

Method:

Column	Ultisil Amy-D, 10µm, 4.6×250mm
Mobile Phase	Methanol
Detection	254nm
Temperature	Ambient
Flow Rate	1.0ml/min
Injection Volume	10ul
Note	/

Chromatogram and Data:

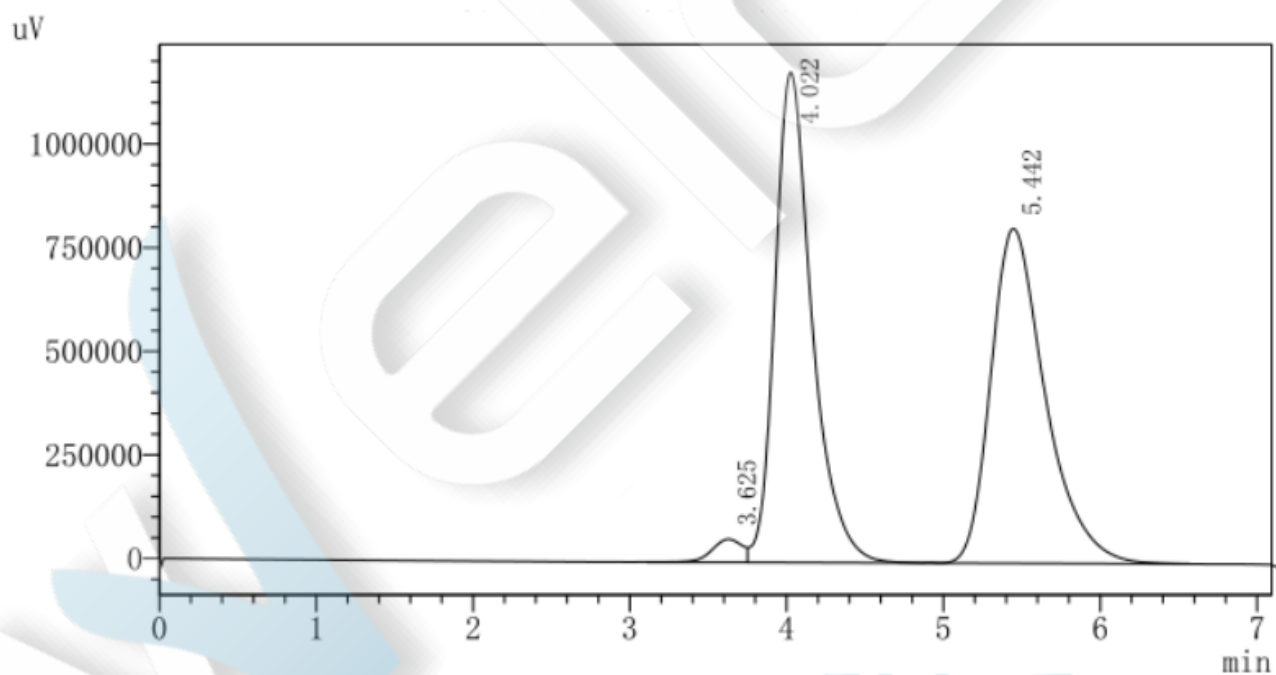

#	RT	Area	Height	Area (%)	Half Width	Plates	Resolution	Tailing
1	3.625	768501	56079	1.93	0.000	1007	0.00	0.000
2	4.022	19527502	1183802	48.98	0.242	1474	0.91	1.373
3	5.442	19568940	807693	49.09	0.361	1228	2.73	1.420

Unless otherwise stated, the results shown in this test report refer only to this sample tested. The report can not be copied without the permission of Welch Materials, Inc.