


## Tromethamine

### Method:

Column	Ultisil® C18-SCX, 5µm, 4.6×250mm P/N: 00212-31043
Mobile Phase	0.025mol/L ammonium dihydrogen phosphate solution
Detection	RID
Temperature	30 °C
Flow Rate	1.0 ml/min
Injection Volume	20 µl
Note	/

### Chromatogram and Data:


#### 1. Midbody

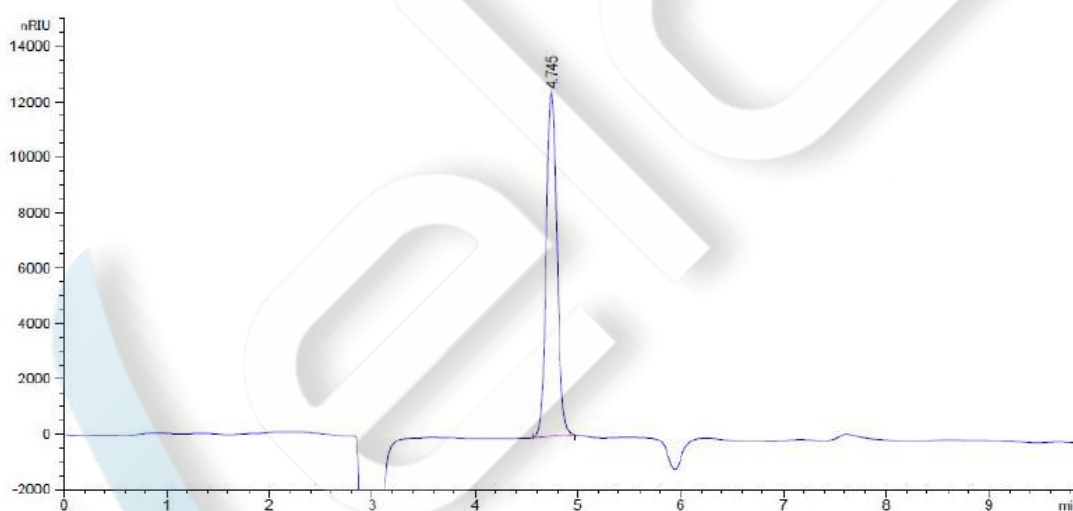

Retention time (min)	K'	Peak Area (mAU*s)	Peak Height (mAU)	Asymmetry	Peak Width (min)	Plates	Resolution	Selectivity
4.745	-	8.83960e4	1.24185e4	0.92	0.1094	10413	-	-


#### 2. Marketable goods

Unless otherwise stated, the results shown in this test report refer only to this sample tested. The report can not be copied without the permission of Welch Materials, Inc.

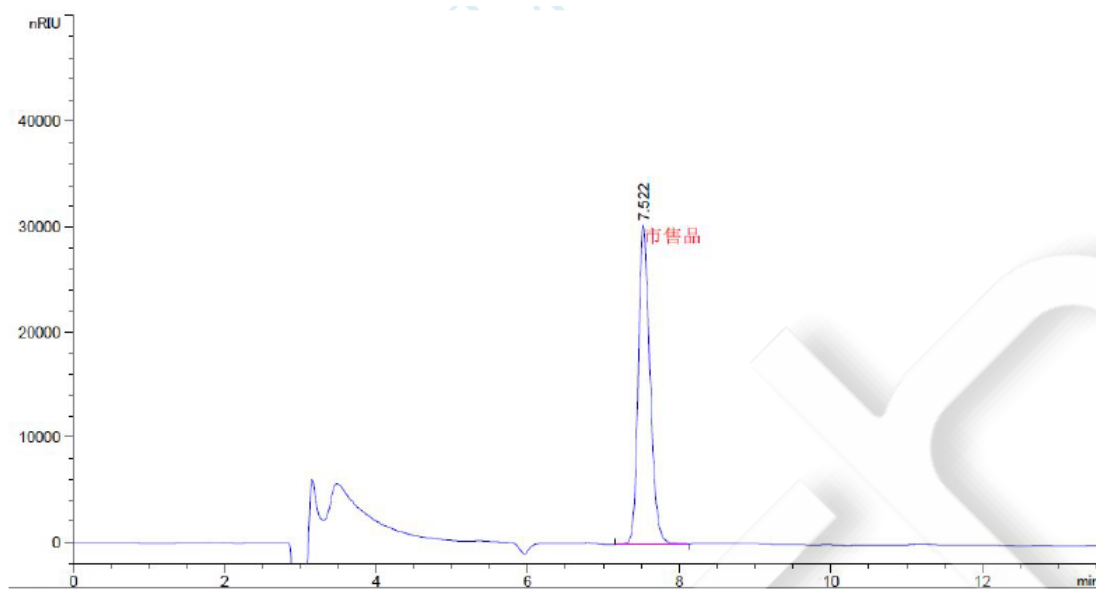

Retention time (min)	K'	Peak Area (mAU*s)	Peak Height (mAU)	Asymmetry	Peak Width (min)	Plates	Resolution	Selectivity
7.522	-	3.30360e5	3.04243e4	0.71	0.1635	11734	-	-

### 3. own product

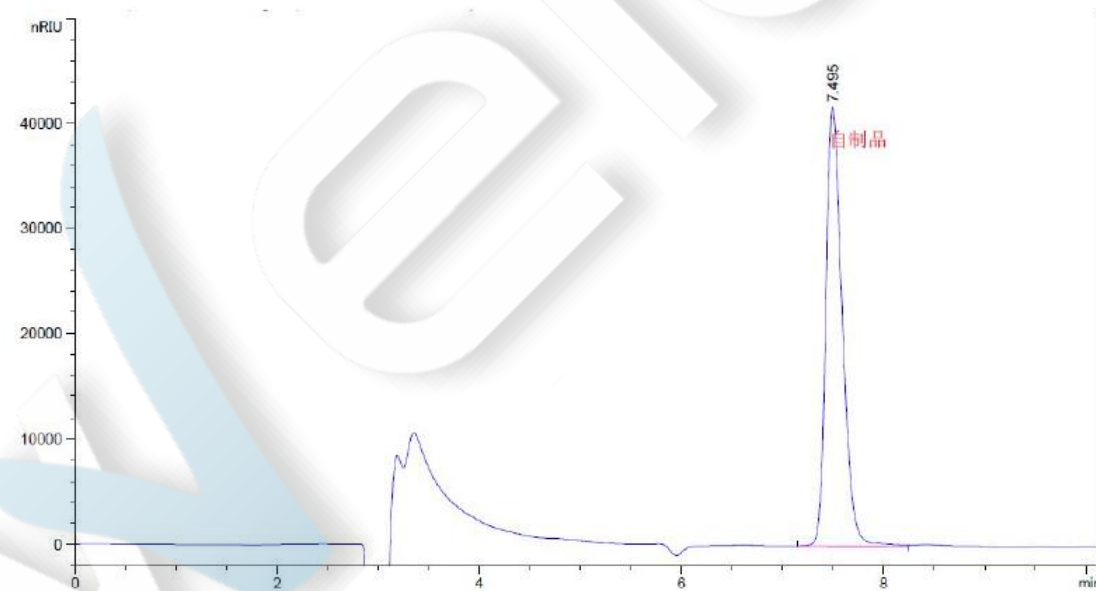

Retention time (min)	K'	Peak Area (mAU*s)	Peak Height (mAU)	Asymmetry	Peak Width (min)	Plates	Resolution	Selectivity
7.495	-	4.73367e5	4.17922e4	0.64	0.1701	10753	-	-

Unless otherwise stated, the results shown in this test report refer only to this sample tested. The report can not be copied without the permission of Welch Materials, Inc.